


Beautifully refined. CROS Pure transmitter.

Based on Signia's proprietary e2e wireless® 3.0, with its industry benchmark for low energy consumption, CROS Pure is the transmitter used for all CROS and BiCROS hearing instrument fittings.



- 1. Directional microphones**
Improve speech understanding in noise by focusing on the person in front of the wearer.
- 2. e2e wireless 3.0**
This third-generation wireless data exchange system enables the energy-efficient audio data transmission to the corresponding *primax* hearing instrument.
- 3. Battery compartment**
For a size 312 battery. With on/off function.
- 4. IP67-rated†**
Resistant to moisture, sweat, dust and dirt.
- 5. Sleek housing**
CROS Pure fits discreetly behind the ear and offers a high degree of wearing comfort.
- 6. miniReceiver 2.0**
Available in four different lengths for retention.
Please note: No output is provided by the receiver when attached to the CROS Pure.

Colors:

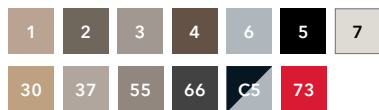

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice.

† Achieved IP67 rating per IEC 60529 standard. Water-resistant device can be completely submerged in water up to 3 feet for 30 minutes and dust will not interfere with the satisfactory operation of the device.

Legal Manufacturer
Sivantos, Inc.
10 Constitution Avenue
Piscataway, NJ 08854
Phone: (800) 766-4500
Fax: (732) 562-6696
www.sivantos.com

Copyright © 2016 Signia GmbH. All rights reserved.
Sivantos, Inc. is a Trademark Licensee of Siemens AG.

signiausa.com SH201662 3.0 2/17 22B771 1K SI/17110-16

signia
Life sounds brilliant.

CROS/BiCROS.

Hear easily - even with unaidable hearing in one ear.



Hearing
Systems

SIEMENS

signiausa.com



Wireless CROS/BiCROS is ideal for patients with an unaidable hearing loss in one ear.

Reduced hearing stress in every situation.

Whether at work in an open-plan office, on the telephone, or at lunch you can hear clearly anytime and anywhere. Thanks to the wireless transmission from the unaidable ear to the ear with better hearing, you can follow conversations effortlessly, even when someone is speaking to you on your unaidable side.

Our new feature, SpeechMaster, identifies the target speech source and isolates it by reducing other competing speech and noise. This innovative technology allows for clear hearing with minimal effort even in a noisy environment.

Enjoy life's true sound.

Signia has just introduced a **CROS/BiCROS** solution for people with an unaidable hearing loss in one ear. Now, not only can you hear and understand clearly, but you can also hear with less strain.

The hearing aids are barely visible and connect wirelessly via the Signia CROS Pure Transmitter. So if someone is speaking to the unaidable side, the signals are processed, transmitted wirelessly and, depending on the level of the hearing impairment in the better ear and the specific environment, the speech signals are amplified.



With Signia's wireless CROS/BiCROS, you can have better hearing everywhere.

Designed for people with unaidable hearing loss in one ear, Signia's CROS and BiCROS hearing solutions offer better hearing all around. A CROS Pure transmitter worn behind the unaidable ear detects sound, processes it, and transmits it to the hearing aid on the other side allowing you to hear sounds from both sides.

